



**Test Fit For Length.** Place on lower teeth only. Keep mouth open.

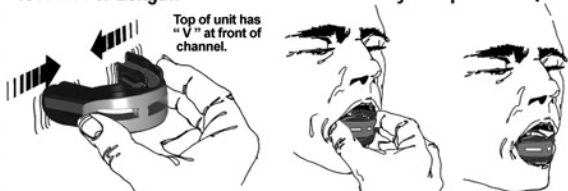


Proper handling for sizing. Lower front teeth determine length.


**Perfect fit.** Lower teeth and molars must fit in channel.

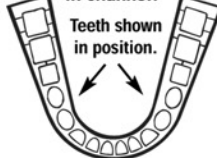

**Slightly too large.** Cut lines (A).

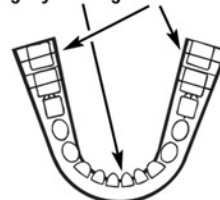

**Too large.** Cut lines (B).

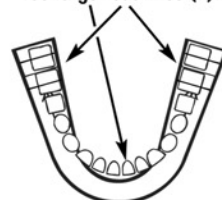


**Trim for Length if needed:** Use larger siccors - Be careful!

Trim/Round all corners for comfort.

**ADULT MODEL:**  
Do not cut in front of line C.

Sizing lines:

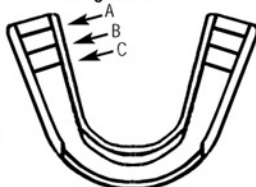

**JUNIOR MODEL:**  
Do not cut in front of line B.




Insert string through the Strap Attachment Hole - located on the side of the post in the breathe space.



**Fully Insert:**  
Breathe space boiling support piece. (included in case).

Carefully insert on the inside of guard.

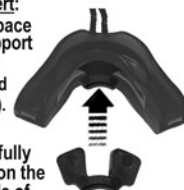

Bring water to rapid boil; then turn off heat to avoid burns.

\* Remove support piece only **AFTER** fitting & cooling guard.



**\* BOILING PROCESS for: PRO+PLUS Models WPR-2004 WPRY-2004 Models Only**

**2 STAGE BOILING PROCESS**

**FIRST STAGE 20 SECONDS**

Place in water up to boil line.

**DO NOT OVER BOIL**



Time critical in seconds. for Models: WPR-2004 / WPRY-2004

**SECOND STAGE 10 SECONDS**

Place in water completely.

**OVER 30 SECONDS TOTAL DESTROYS THE WPR-2004**




**\* BOILING PROCESS for: PRO+PLUS HIGH PERFORMANCE WPR-300 Model Only**

**2 STAGE BOILING PROCESS**

**FIRST STAGE 30 SECONDS**

Place in water up to boil line.

**DO NOT OVER BOIL**

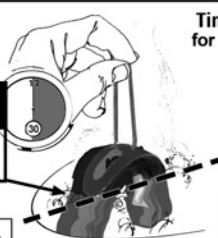

Time critical in seconds. for HIGH PERFORMANCE Model: WPR-300

**SECOND STAGE 15 SECONDS**

Place in water completely.

**OVER 45 SECONDS TOTAL DESTROYS THE HP WPR-300**




Shake off excess water.

Leave Strap Attachment inserted.

Lower front flange may curl in.

Use thumbs to open curled channels (check Top and Bottom) before Molding.

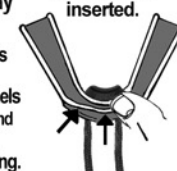

Leave Strap Attachment and insert in place to mold!

Open mouth, slightly extending the lower jaw forward.

Place on lower teeth.




Bite down slowly to insure teeth fit in upper channel.

Bite down to imprint teeth.

**! CAUTION !**

IF you have BRACES, BRACKETS or APPLIANCES on your teeth, **SKIP** the following BOXED Step.

**DO NOT MOLD OUTSIDE SURFACES INTO BRACES BRACKETS or APPLIANCES !**

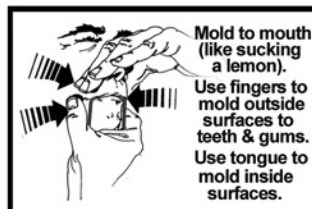

Remove & rinse with cool water.

Remove breathe space boiling support piece.

**Check fit.**




Carefully Re-boil if necessary.

**PROPER STRAP ATTACHMENT: IMPORTANT !:**

**DO NOT TIE or WRAP ANY part of attachment strap or string !** Doing so, can negate Safety Release clasp's effectiveness.



**Optional:**  
Attach 2nd segment (included in case) to increase length for Lanyard style use.

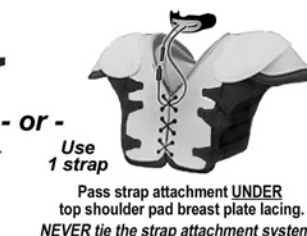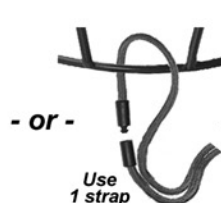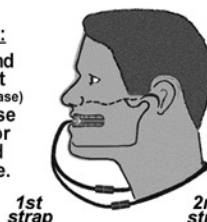

Pass strap attachment **UNDER** top shoulder pad breast plate lacing. **NEVER** tie the strap attachment system !

## CARE & CLEANING:

- Keep your Brain Pad & storage case clean. Store the Brain Pad in its case when not in use.
- Wash the Brain Pad after each use in lukewarm water. Hot water may soften & distort the Brain Pad.
- Sanitize regularly with an astringent mouthwash for good hygiene.
- With proper care, your Brain-Pad should provide protection for a whole season.